



User Manual

Breezy Tag (SKU: BTAG01)



Manufacturer and Contact Information

Manufacturer: Breezy Med

*Address: 317 Commercial St NE, Suite 101
Albuquerque, NM 87102*

For product information, sales, service, and technical support:

<https://breezymed.com/>

Customer Support Phone: (833) 273-3991

Email: info@breezymed.com

This manual covers the model: Breezy Tag (BTAG01)

SPECIFICATIONS:

Wireless protocol	Bluetooth® BLE 5.0
Size	Diameter: 1.02 in. (26.0 mm), Width: 0.28 in (7.1 mm)
Weight	0.2 oz. (5 g)
Certifications	FCC CE RoHS REACH
Battery	CR2032 220mAh Replaceable
Battery Life	Approximately 1.9 years
Mounting method	3M VHB tape
Operating temperature range	-4°F to 140°F
Operating humidity range	0 to 95% relative humidity, non-condensing
Breezy Tag location - operating distance	Within 18.5 feet of Breezy Smart Controller during fogging
Required equipment for asset tracking	Breezy Smart Controller (Gen 2), Breezy Blue, Breezy Buggy


Environmental Considerations

This section provides information about the environmental impact of the product.

Product End-of-Life Handling

Observe the following guidelines when recycling an instrument or component:

Equipment recycling. *Production of this device required the extraction and use of natural resources. The device may contain substances that could be harmful to the environment or human health if improperly handled at the product's end of life. To avoid the release of such substances into the environment and to reduce the use of natural resources, we encourage you to recycle this product in an appropriate system that will ensure that most of the materials are reused or recycled appropriately. It is recommended that this device be taken to a certified Waste Electronic and Electrical Equipment recycler.*

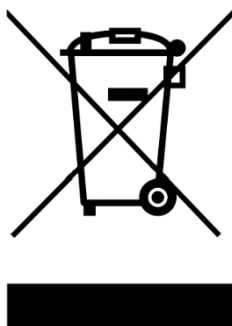

Do not throw in trash at end of product life.

Contents


Manufacturer and Contact Information	2
Environmental Considerations	3
Breezy Tag Setup.....	5
Installing Breezy Tag	6
Replacing the Battery	7
Breezy Tag management.....	8
Troubleshooting	9
LIMITED WARRANTY (USA/CANADA).....	10

Breezy Tag Setup

1. Remove the plastic battery insulating pull tab from the side of the Breezy Tag. A red light will flash quickly on the side of the Breezy Tag displaying the QR Code.



2. After approximately 5 seconds, the Breezy Tag should be advertising
3. Using a Breezy Smart Controller that has been configured to the site, tap the Breezy Tag button on the devices Screen.

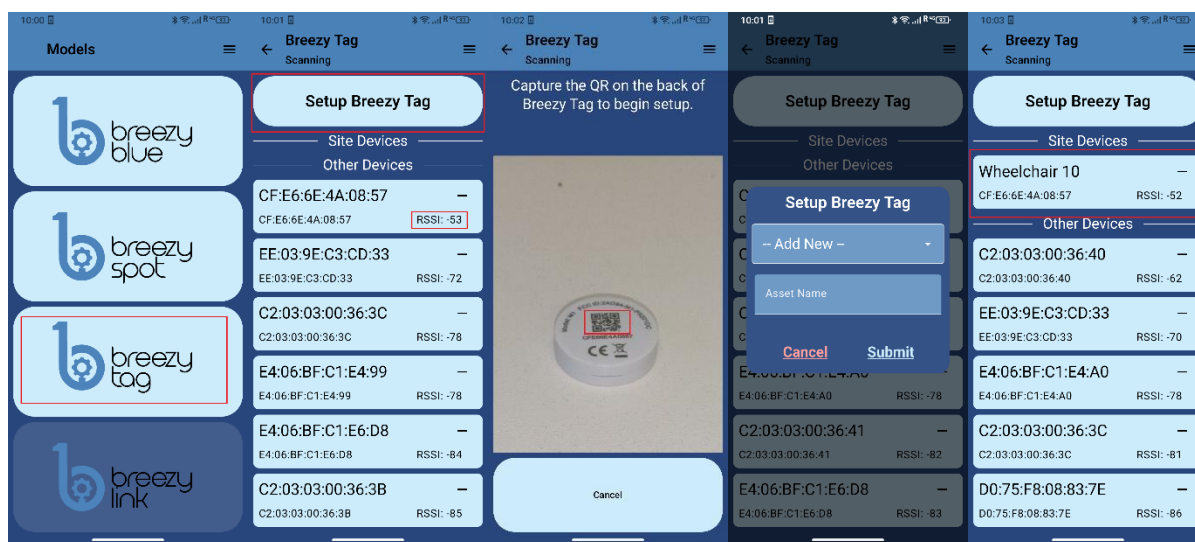

4. When the controller is stationary, the list will be populated with Breezy Tag devices within Bluetooth range. If the controller is in contact with a Breezy Tag, that Breezy Tag will report a strong RSSI value (generally -30 to -55).

5. Tap the Setup New Breezy Tag button.
6. Using the camera on the rear of the controller, move the QR code on the bottom of the Breezy Tag into the image shown on the screen. When the QR code has been read, the camera will close. Move the controller closer and further from the QR code if you have issues reading the QR code.
7. Enter a name for the asset (example: Wheelchair 10) and press continue.
8. The Breezy Tag device will be shown in the Site devices and named with the asset name,
9. Asset names and Breezy Tags associated with an asset can be updated at <https://cloud.breezymed.com/assets>

Installing Breezy Tag

1. Select an asset to mount the Breezy Tag device. It should be near the location of the Breezy Blue during disinfection.
 - a. Breezy Tag should be within 18.5 feet of the Breezy Smart Controller during disinfection.
2. Ensure the mounting location is free of dust or other contaminants that may prohibit adhesion.
3. Remove the red adhesive liner from the double stick 3M tape and firmly press against Breezy Tag and the selected asset for at least 5 seconds.

Replacing the Battery

1. Breezy Tag battery levels can be monitored 3 ways
 - a. The battery state can be found at <https://cloud.breezymed.com/devices> and selecting Breezy Tag
 - b. Breezy Tag low battery notifications can be sent to email or text by enabling notifications at <https://cloud.breezymed.com/settings>
 - c. Breezy Tag will flash a red indicator light twice every 10 seconds when battery is low
2. The battery can be replaced by
 - a. Use a small flat blade screwdriver to pry the top off the Breezy Tag.
 - b. Remove the CR2032 coin cell battery and replace it.
 - c. Snap the top back on by aligning the depression on the side of the top plate with the slot on the bottom of the bottom plate, then pressing firmly straight down.



- d. Once the battery is replaced, a red light will quickly flash.

Breezy Tag management

1. Breezy Tag and the location name will be shown on Breezy Cloud at <https://cloud.breezymed.com/devices> and the asset will be recorded in the Breezy Blue fog log when disinfecting.
2. Breezy Tags can be removed from the site at <https://cloud.breezymed.com/devices>

Requirements for Successful Asset Disinfection

1. For a Breezy tag to report an asset as successfully disinfected, it must meet these criteria:
 - a. Breezy Blue completed a full fog event was allowed its full time to settle.
 - b. The Smart Controller was present in the room and stationary during the fogging as settling events
 - c. The asset with the Breezy Tag attached was present in the room and stationary during the fogging and settling events.
2. If all the criteria are met, the asset will be considered disinfected on the Cloud site

Assets								Advanced
Asset Name	Fog Ok ⓘ	Settling Time ⓘ	Controller Stationary ⓘ	In-Room ⓘ	Present Continuously ⓘ	Stationary ⓘ	Asset Disinfected	
wheelchair 2	✓	✓10 minutes	✓	✓	✓	✓	✓	
wheelchair 1	✓	✓10 minutes	✓	✓	✓	✓	✓	
insulin pump 1	✓	✓10 minutes	✓	✓	✓	✓	✓	

Troubleshooting

PROBLEM	POSSIBLE CAUSES	SOLUTION
<i>Unable to add Breezy Tag to site using Breezy Smart Controller</i>	<i>Breezy Tag is not powered ON</i>	<i>Make sure the plastic pull tab has been removed.</i>
<i>Asset is not marked as disinfected on Breezy Cloud</i>	<i>Not all the criteria were met when the disinfection even occurred</i>	<i>Follow Steps 1-2 under "Requirements for Successful Asset Disinfection"</i>
<i>Assets aren't shown in the fog log even though a Breezy Tag is present</i>	<i>Breezy Smart Controller was moving or too far away from Breezy Tags during fogging.</i>	<i>Ensure Breezy Tags and the Breezy Smart Controller are in proximity and aren't moving during fogging</i>

LIMITED WARRANTY (USA/CANADA)

Breezy Med warrants that this product will be free from defects in materials and workmanship for a period of one (1) year from the date of shipment. If any such product proves defective during this warranty period, Breezy Med, at its option, either will repair the defective product without charge for parts and labor, or will provide a replacement in exchange for the defective product. Parts, modules and replacement products used by Breezy Med for warranty work may be new or reconditioned to like new performance. All replaced parts, modules and products become the property of Breezy Med.

In order to obtain service under this warranty, Customer must notify Breezy Med of the defect before the expiration of the warranty period and make suitable arrangements for the performance of service. Customer shall be responsible for packaging and shipping the defective product to the service center designated by Breezy Med, with shipping charges prepaid. Breezy Med shall pay for the return of the product to Customer if the shipment is to a location within the country in which the Breezy Med service center is located. Customer shall be responsible for paying all shipping charges, duties, taxes, and any other charges for products returned to any other locations.

This warranty shall not apply to any defect, failure or damage caused by improper use or improper or inadequate maintenance and care. Breezy Med shall not be obligated to furnish service under this warranty a) to repair damage resulting from attempts by personnel other than Breezy Med representatives to install, repair or service the product; b) to repair damage resulting from improper use or connection to incompatible equipment; c) to repair any damage or malfunction caused by the use of non-Breezy Med supplies; or d) to service a product that has been modified or integrated with other products.

THIS WARRANTY IS GIVEN BY BREEZY MED WITH RESPECT TO THE PRODUCT IN LIEU OF ANY OTHER WARRANTIES, EXPRESS OR IMPLIED. BREEZY MED AND ITS VENDORS DISCLAIM ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. BREEZY MED' RESPONSIBILITY TO REPAIR OR REPLACE DEFECTIVE PRODUCTS IS THE SOLE AND EXCLUSIVE REMEDY PROVIDED TO THE CUSTOMER FOR BREACH OF THIS WARRANTY. BREEZY MED AND ITS VENDORS WILL NOT BE LIABLE FOR ANY INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES IRRESPECTIVE OF WHETHER BREEZY MED OR THE VENDOR HAS ADVANCE NOTICE OF THE POSSIBILITY OF SUCH DAMAGES.